www.coretechortho.com

coretech® by vivehealth

OWNER'S MANUAL
360 EXO WALKER TALL

SKU: SUP2032BLK



To see all FAQs in one place, visit, **coretechortho.link/0ha**


SUP2032BLK

360 EXO WALKER TALL

PUTTING THE WALKER ON

For a video demonstration
check out **coretechortho.link/4h6**

1. Unfasten all straps and open the liner.
2. (Optional) Two extra pads are included to cushion the inside of the liner, they can be used around the ankle, the calf, under the foot, or the straps.

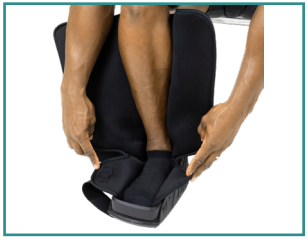



3. With your foot inside the liner secure the bottom strap, then the top. Make sure your heel is positioned completely back.





4. Secure the outer straps from bottom to top and tighten to a comfortable fit. Do NOT overtighten. Make sure the outer plastic splint is centered along the front arch of the foot and the velcro tabs are secured at the top.





5. Set the air valve to “INFLATE” by turning it clockwise.



6. Press the pump to inflate to your desired compression.



7. Set the air valve to “DEFLATE” by turning counter-clockwise, use hand compression to speed up the deflating process.



8. To remove or adjust the toe cover, unpin the hooks on each side by positioning them vertically and pushing in to unlock. Remove the cover completely or reposition for a better fit.





CARE INSTRUCTIONS

- Remove the padded liner.
- Spot clean with a damp cloth. Do NOT submerge in water. Do NOT machine wash or dry. Air dry.
- Wipe the plastic with a damp cloth and dry immediately.

www.coretechortho.com

coretech®

by vivehealth

Distributed by

vivehealth

8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com

261223