www.coretechortho.com

coretech®
by vivehealth

OWNER'S MANUAL
ULTRA 639 LSO

SKU: SUP3114BLK




To see all FAQs in one place, **visit coretechortho.link/iak**

SUP3114BLK

ULTRA 639 LSO

WHAT'S INCLUDED


- 1x 639 LSO Ultra
- 2x Lateral Side Panels (Removable)
- 2x Extenders
- 1x Adjustable Anterior Panel

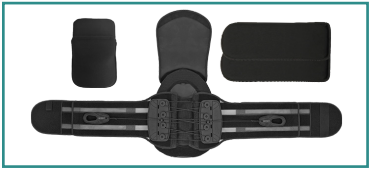


SIZING GUIDELINES

Before application, ensure proper sizing by:

- Measuring waist circumference around the belly button.
- Referencing the size chart printed inside the belt.
- Adjusting lateral panels as needed for optimal fit.
- Extending or retracting the anterior panel (up or down) for appropriate stabilization.
- Testing compression levels before securing pull tabs.

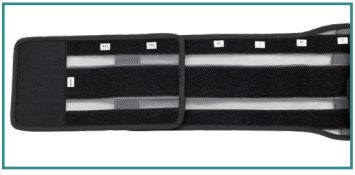

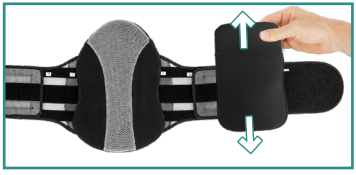

ADJUSTING THE BRACE

1. Disconnect the side panels and size to the waist needed. Use the included extenders if necessary. Size each side evenly.
2. Lateral panels can be secured to each side of the belt for support if required.





3. The anterior panel uses hook and loop to secure to the inside of the belt on the left hand side. It can be repositioned as needed.



APPLICATION

Assistance may be required for proper application!

1. Position the back panel centered on your spine so that it covers the S1 up to T9 vertebrae. It should cover the entire lower back.




2. Position the anterior panel and center it on the abdomen. Assistance may be required. Overlap the left side belt over and then secure the right side to complete the fit.





3. Pull the rubber pulley tabs to compression, do not over tighten.

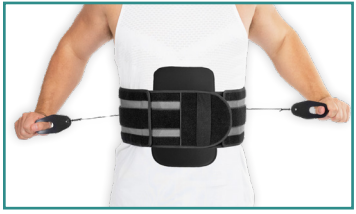

The anterior panel should be positioned at the center of your abdomen with the curved side at the bottom and the straight edge at the top.



Test fit and ensure compression is not too tight or too loose, adjust as needed.



CARE INSTRUCTIONS

Separate cover from plastic splint and pulley system. Hand wash and let air dry. Do not machine wash or tumble dry.

WARNINGS ⚠

This product is intended for application by a licensed healthcare practitioner as directed by a physician or other qualified medical authority. Permanent injury may result from use other than as directed by a physician or other medical authority.

TROUBLESHOOTING

- **If the brace feels too loose:** Increase tension on the pull tab cords or adjust to a smaller size setting.
- **If the brace feels too tight:** Reduce tension on the pull tab cords or adjust to a larger size setting.
- **If the front panel doesn't align properly:** Ensure the back panel is centered on the spine and readjust the lateral panels.

WARRANTY

The 639 LSO Ultra is covered by a 1-year warranty against all mechanical or manufacturing defects, from the date of purchase to the original owner.

For additional support and frequently asked questions, visit: **coretechortho.link/iak**

NOTES

www.coretechortho.com

coretech®

by vivehealth

Distributed by

vivehealth

8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com
Made in China

230626